
DELUXE MASSAGE CHAIR

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Please read the user manual carefully before use.

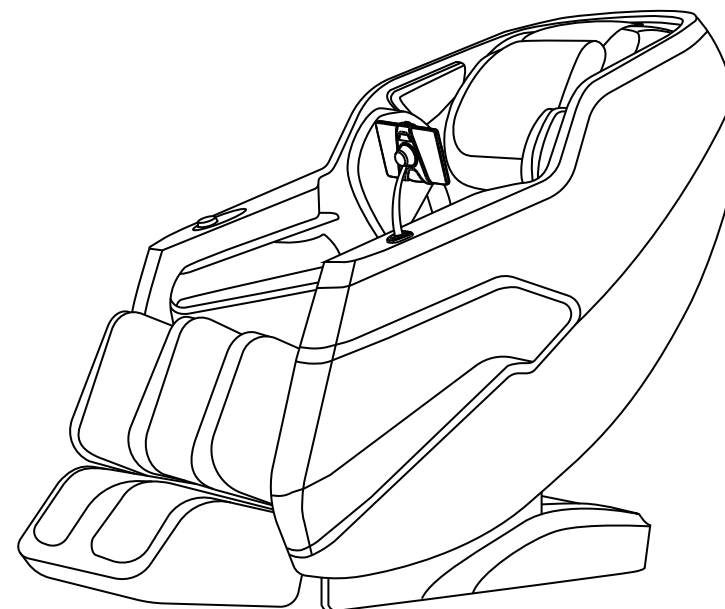Thank you for choosing our product and your support. In order to understand the directions fully, please read user's manual carefully. And after reading, please keep it for next browsing. We reserve the right of final explanation to the design of product. Dimensions and specifications are subject to change without notice for further improvement.

Contents

Safety	1-6
Preparation Before Use	7-11
Move the Chair	7
Site of Placement	7
Install Legrest	8
Install Backrest Pad and Pillow Pad	9
Installation of Manual Bracket	10
Remote Controller Connection	10
Connect Power	11
Bluetooth pairing	11
Components	12-14
Outer Structure	12
Controller	13
Shortcut Key	14
Corpus Summary	15-16
Cleaning and Maintenance	17
Troubleshoot	18
Specifications	19

Safety

IMPORTANT SAFETY INSTRUCTIONS

This massage chair is intended for household only .When using an electrical appliance, basic precautions should always be followed, including the following:

Read all instructions before using this massage chair.

DANGER-To reduce the risk of electric shock:

- 1) Always unplug this appliance from the electrical outlet immediately after using and before cleaning.
- 2) Never use pins or other metallic fasteners with this appliance.
- 3) Carefully examine the covering before each use. Discard the appliance if the covering shows any sign of deterioration, such as checking, blistering, or cracking.
- 4) Keep Dry - Do not operate in a wet or moist condition.

WARNING-To reduce the risk of burns, fire, electric shock, or injury to persons:

- 1) An appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
- 2) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- 3) Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons. Keep children away from extended foot support.
- 4) Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- 5) Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water.
Return the appliance to a service center for examination and repair.
- 6) Do not carry this appliance by supply cord or use cord as a handle.
- 7) Keep the cord away from heated surfaces.
- 8) Never operate the appliance with the air openings blocked. Keep the air openings free of lint, hair, and the like.
- 9) Never drop or insert any object into any opening.

Safety

- 10) Do not use outdoors.
- 11) Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- 12) To disconnect, turn all controls to the off position, then remove plug from outlet.
- 13) Temperatures sufficiently high to cause burns may occur regardless of the control setting. Do not use on an infant or invalid or on a sleeping or unconscious person. Do not use on insensitive skin or on a person with poor blood circulation. Check the skin in contact with the heated area of the appliance frequently to reduce the risk of blistering.
- 14) Do Not Crush - Avoid sharp folds.
- 15) Do not use massager in close proximity to loose clothing or jewelry.
- 16) Keep long hair away from massager while in use.
- 17) Connect this appliance to a properly grounded outlet only. See Grounding Instructions.
- 18) The apparatus shall not be exposed to dripping or splashing and that no objects filled with liquids, such as vases, shall be placed on the apparatus.
- 19) An apparatus with CLASS I construction shall be connected to a MAINS socket outlet with a protective earthing connection
- 20) The appliance coupler is used as the disconnect device, the disconnect device shall remain readily operable.
- 21) This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- 22) If it is not in use for a long time or after use, please turn off the power and pull out the power plug.
- 23) Do not damage, break, over-bend, compress, tighten, wind, rotate or knot the power cord.

Safety

SAVE THESE INSTRUCTIONS

This massage chair is intended for household use only.

● Usage Environment

- Do not use the massage chair in an excessively humid or dusty environment as this may result in malfunction or electric shock.
- Do not use the massage chair in a room with a temperature of 104°F or higher.
- Do not expose the massage chair to heaters, stoves or direct sunlight.
- Place and use the massage chair on an even, non-slip floor.
- The massage chair is designed for home use and not for commercial purposes.

Caution

- The massage chair is to be used indoors only.
- If you use the massage chair in a cold room, do not increase the room temperature abruptly. It is recommended that you increase the temperature gradually to a normal level.
- If the massage chair is stored in a cold place and is being brought into a warm environment, it is advisable to wait one hour before use. This is because its performance may not be at optimum level, as there may be water droplets on the mechanical parts due to condensation. Using the massage chair under such conditions may result in malfunction.

IMPORTANT SAFETY INSTRUCTIONS

● Safety Precautions

- To ensure safe and correct use of the massage chair, do not operate it without reading this operation manual carefully.
- Do not use the massage chair when the upholstery or seat is removed. This may cause injury or malfunction.
- Do not use the massage chair in combination with other therapeutic equipment or electric blanket, etc., as this may result in ineffectiveness or injury.
- Do not use the massage chair an hour before or after eating as this may have adverse effects on the user.
- Do not allow children or pets to play around the massage chair, e.g., behind the backrest, under the seat or leg rest of the massage chair. This may result in injury.

Safety

- Do not rest or place heavy objects on the armrests, leg rest or backrest as this may result in malfunction or injury.
- Do not use the massage chair when your body is wet.
- Do not operate the massage chair with wet hands.
It is not recommended to use the massage for more than 30 minutes.
- Do not perform continuous massage on the same spot of your body for more than 5 minutes at a time as this may result in excessive stimulation and may have adverse effects.
- Do not insert your hand or foot along the paths of the massage rollers during use as this may cause injury.
- Stop operation immediately if you experience any discomfort.
- This product is not intended for self-treatment of conditions that should be managed by a qualified health care provider.
- People on medication or with medical conditions, please consult your doctor before use.
- Please do not use the product if you are not feeling well.
- Remove all scarves, neckties, necklaces and jewelry before using the massager.

● People with Medical Conditions

Consult your doctor before use if you:

- are under medical rest as ordered by a doctor.
- have spinal disorders, an abnormal spinal condition or have suffered a spinal injury.
- have back problems.
- have diabetes, osteoporosis or sensory impairment.
- have joint dysfunction.
- have a pacemaker or other electronic medical devices.
- may be pregnant.
- have phlebitis or thrombosis.
- have an increased risk for blood clots.
- had recent surgery.
- have surgical pins, screws, or anything mechanical implanted in your body.

Safety

● Instances When the massage chair Should Not Be Used

In the following events, please switch off the main power switch immediately and disconnect the power cord from the electrical outlet:

- If water is accidentally spilled onto the massage chair. This may result in electric shock or malfunction.
- If the fabric of the backrest is torn and the internal components are exposed.
- If pain or discomfort is felt during massage, stop operation immediately and consult your doctor.
- If you detect any malfunction or any other abnormal condition during operation.
- If there is a power failure. Injury may occur when power is restored unexpectedly.
- If there is lightning.

● Assembly and Repair of the massage chair

- Do not disassemble the backrest cover of the massage chair. Touching the internal components may result in malfunction or electric shock.
- The massage chair must only be professionally serviced or repaired. Do not attempt to disassemble or repair the massage chair yourself.

● Things to Note About the Power Plug and Cord

- Check that the Alternating Current (AC) voltage corresponds with the specifications indicated for the massage chair.
- Do not connect or disconnect the power plug from the electrical outlet with wet hands. This may result in malfunction or electric shock.
- When disconnecting the power plug, pull it by holding the plug, not the cord.
- Although operation will be stopped automatically by the Auto-Timer function, always remember to switch off the main power switch after use.
- Do not place the power cord under the massage chair or any other heavy objects.
- Do not wind the power cord around the massage chair as this may damage the cord and result in a fire or electric shock.
- Do not operate the massage chair with a damaged cord or extension cord.
- Do not use the massage chair if the electrical outlet is loose.
- If the power cord is damaged, it should be replaced by a qualified technician.


Safety

● Grounding Instructions

This product must be grounded. If it should malfunction or breakdown, grounding provides a path of least resistance for electric current to reduce the risk of electric shock. This product is equipped with a cord having an equipment-grounding conductor and a grounding plug. The plug must be plugged into an appropriate outlet that is properly installed and grounded in accordance with all local codes and ordinances.

DANGER - Improper connection of the equipment-grounding conductor can result in a risk of electric shock. Check with a qualified electrician or serviceman if you are in doubt as to whether the product is properly grounded. Do not modify the plug provided with the product - if it will not fit the outlet, have a proper outlet installed by a qualified electrician.

This product is for use on a nominal 100-240V~ circuit and has a grounding plug that looks like the plug illustrated in sketch A as following illustration. Make sure that the product is connected to an outlet having the same configuration as the plug. No adapter should be used with this product.

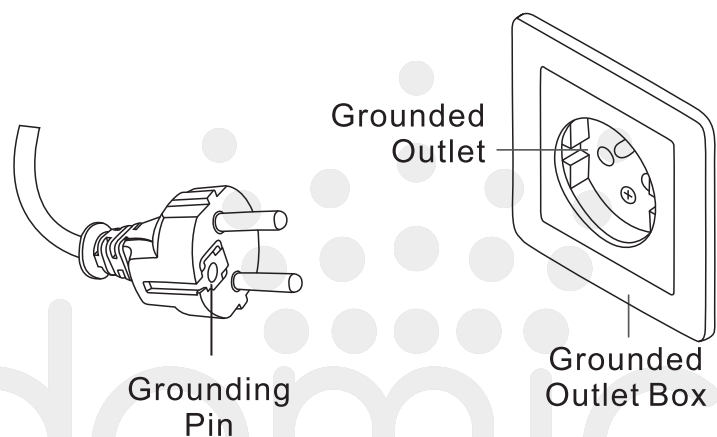
06

Preparation Before Use


Move the Chair

- Please turn off the main power and unplug the power from the electrical outlet.
- Tilt the chair backward as shown in the picture below and wheel it to the desired location.
- Please cover the floor with mat when move the massage chair in case of any damage of floor will be made.

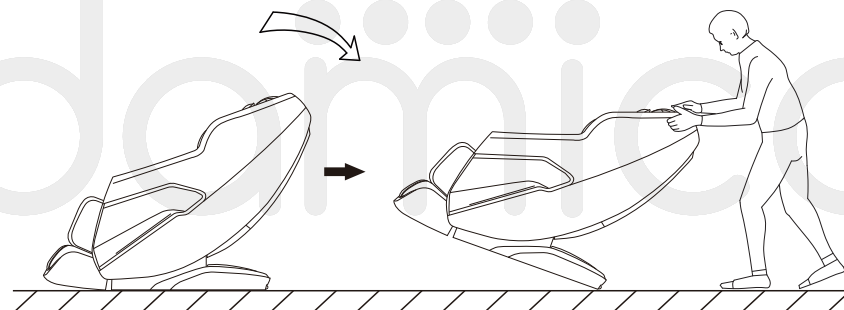Site of Placement

Please make sure that there is an allowance of more than 5cm behind the massage chair, and more than 35cm in front of the massage chair to allow the backrest to recline without obstruction during operation.


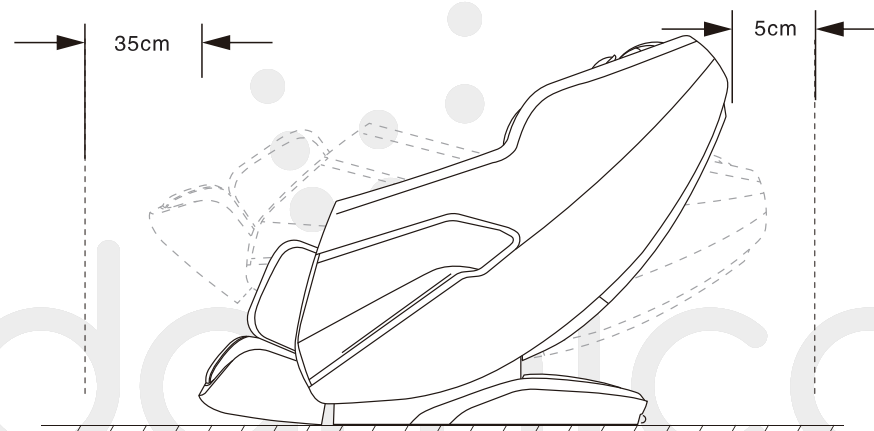
07



Preparation Before Use

Install Legrest




- Attach the push rod of the legrest to the frame, then insert the pin to make it locked.
- Insert the air pipe ① and interface ② into the counterpart position of the chair, separately.

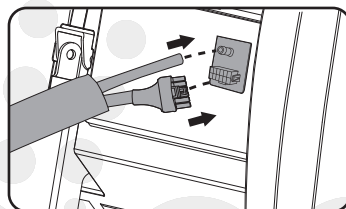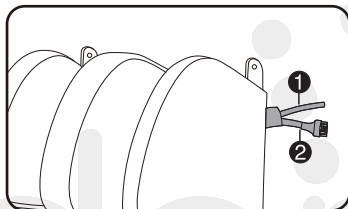- Put the legrest bracket ③ into the fixed axis ④.

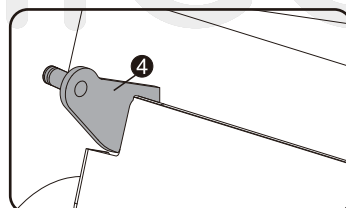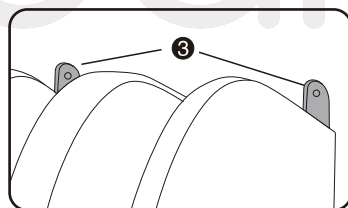- Lock the footrest by the clips, screw it tightly.

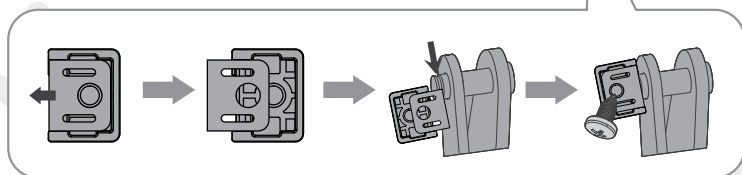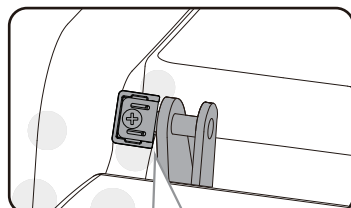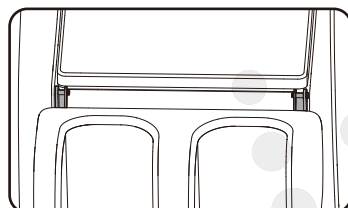08


Preparation Before Use

Install Backrest Pad and Pillow Pad

- Attach the upper backrest pad to backrest with zipper.
- Attach the pillow pad to backrest with velcro.

- Attach the left side of backrest pad to backrest with zipper, same as the other side.

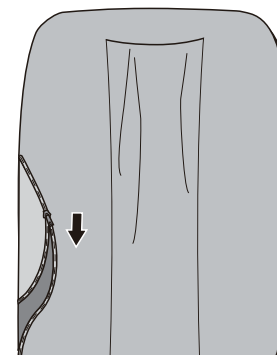
09

Preparation Before Use

Installation of manual bracket


- Aligned holes, lock screw.
- Press down to align the snap.

Note: The arrow direction is towards the front of the handrail.

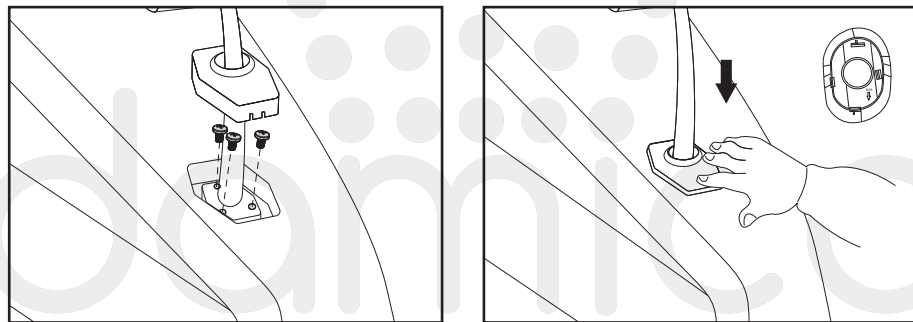Remote Controller Connection

- Connect the cable of remote controller with the cable under seat.

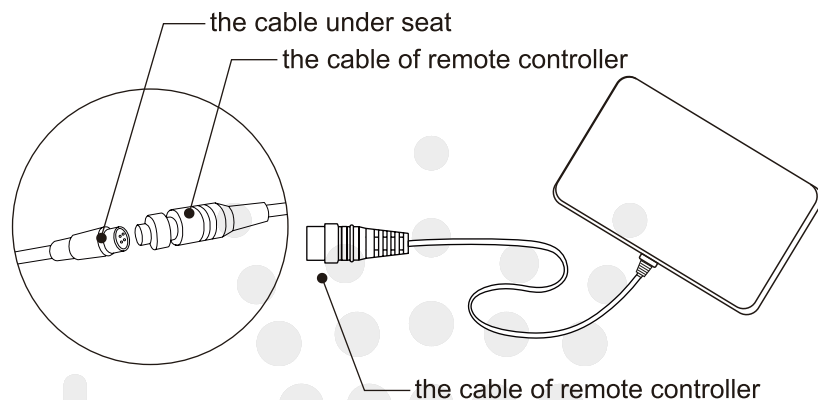
10

Preparation Before Use


Connect power

- Plug the power cord into power inlet, turn on the power switch.

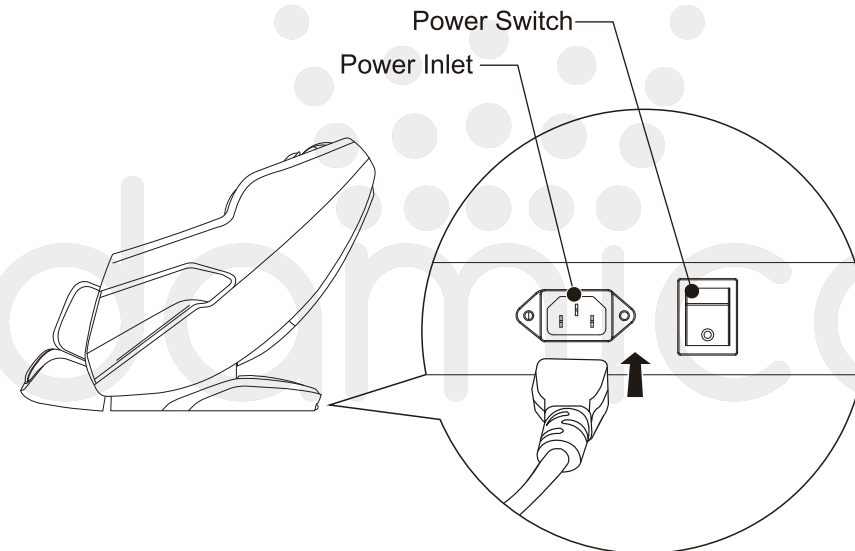Bluetooth pairing



On your phone Bluetooth list, choose the model "Damico Altro 3D" to do the pairing.

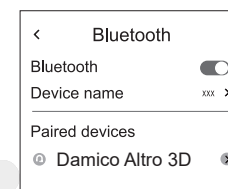11

Components

Outer Structure

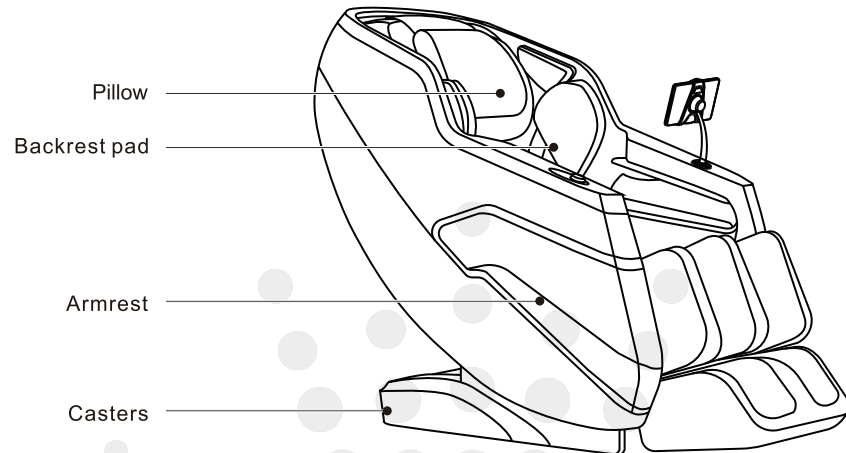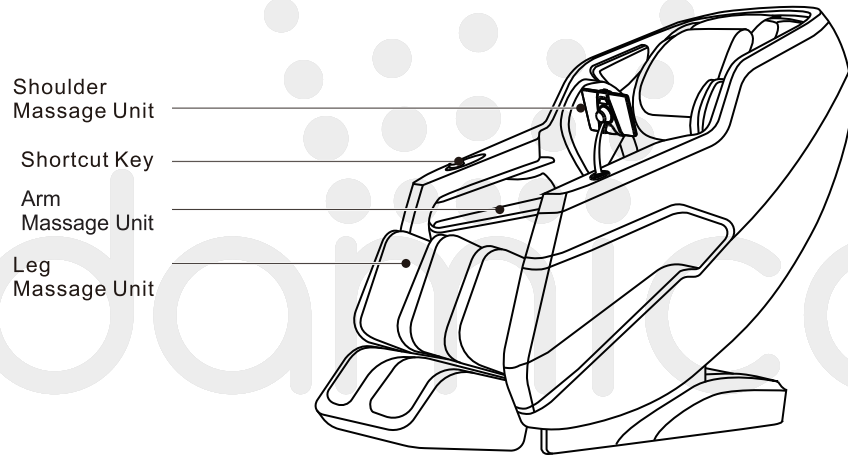12




Components

Controller

(The specific display is subject to the actual product.)



Shoulder & Neck Testing

- Press and hold the power button on the side of the hand controller to turn on the massage chair.
- After powering on, select the massage program and initiate body shape detection. Please relax and lie down comfortably.
- If the body shape detection is not in the optimal position, press the () () button and then press "OK" to start the massage.

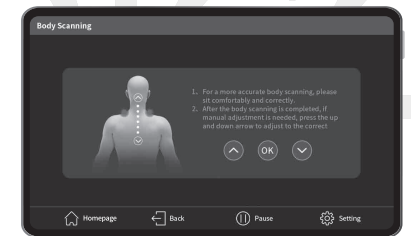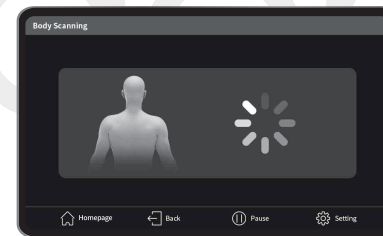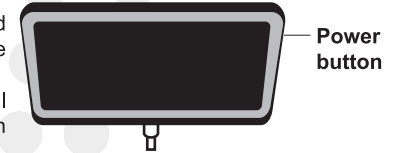Function Introduction

- You can select the desired functions based on your personal needs.

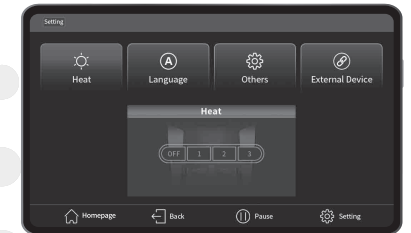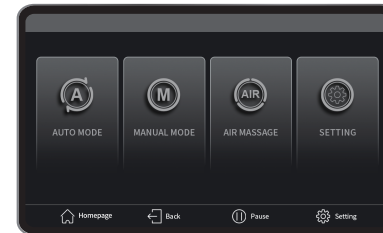Auto Program: Enter the automatic program interface to choose from a variety of massage programs.

Manual Program: Enter the manual program interface to choose massage techniques, massage areas.


Airbag Massage: Enter the airbag massage interface to select airbag massage areas and zero gravity levels.

Setting: Enter the setting interface to configure features such as heating, language, and Bluetooth.

13

Components

Shortcut Key

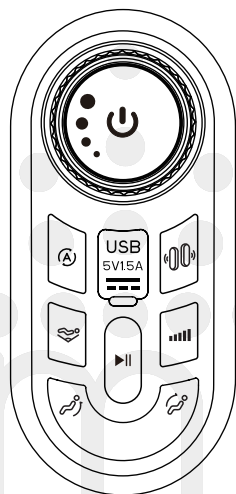Armrest shortcut controls

FUNCTION	ICON	DESCRIPTION
Power		Press this button to activate or deactivate the massage chair.
Rotation		Rotate clockwise to increase stage, and counterclockwise to decrease stages. The stage are divided into 1-6 stage.
Auto		Press this button to switch between multiple automatic program functions. The specific display is subject to the actual product.
Pause		Press this button to pause all massage functions. Press again to resume the massage.
Airbag		Press this button to turn on / turn off airbag massage function / switch airbag position.
Intensity		Press this button to adjust the strength of the airbag, cycling in order from small to large.
Zero-G		Press repeatedly this button to activate one of the three zerogravity modes.
Back Up		Long press this button to adjust the up of the backrest, and release the button to lock the current angle of the backrest.
Back Down		Long press this button to adjust the down of the backrest, and release the button to lock the current angle of the backrest.
USB	USB	You can remove the protective cover (tab) to recharge your phone using a USB cable (not provided).

Corpus Summary

Corpus Summary

Intention	Corpus	Reply
Basic Functions	Hello massage chair	Hello Master
	Power on	Please select your favorite program to start your massage.
	Power off	During shutdown, please leave or restart the massage chair after the chair is completely reset. Thank you and hope you enjoy the massage.
	Pause massage	Massage paused
	Continue to massage	Massage resumed
	Turn on Bluetooth	OK, Bluetooth on
	Turn off Bluetooth	OK, Bluetooth off
	Turn off Foot roller	OK, Foot roller off
	Turn on Foot roller	OK, Foot roller on
	Turn on Heating	OK, Heating on
	Turn off Heating	OK, Heating off
	Zero gravity	OK,Command complete
	Back up	OK,Command complete
	Back down	OK,Command complete
	Leg up	OK,Command complete
	Leg down	OK,Command complete
Start airbag massage	OK, Start airbag massage	
Stop airbag massage	OK, Stop airbag massage	
Auto Massage setting	Master massage	Start Master massage
	Relaxing massage	Start Relaxing massage
	Deep Shiatsu massage	Start Deep Shiatsu massage
	Refreshing massage	Start Refreshing massage
	Wake Up massage	Start Wake Up massage
	Noon Rest massage	Start Noon Rest massage
	Neck & Shoulder massage	Start Neck & Shoulder massage
	Waist & Back massage	Start Waist & Back massage
Slim massage	Start Slim massage	

Corpus Summary

Intention	Corpus	Reply
Auto Massage setting	Sole Relieving massage	Start Sole Relieving massage
	Spine Care massage	Start Spine Care massage
	Calf Care massage	Start Calf Care massage
Miscellaneous Function	Volume up	OK, Volume increase
	Volume down	OK, Volume decrease
	Maximum volume	OK, Command complete
	Minimum volume	OK, Command complete
	Kneading massage	Kneading massage
	Shiatsu massage	Shiatsu massage
	Knocking massage	Knocking massage
	Scrapping massage	Scrapping massage
	Tapping massage	Tapping massage
Clapping massage	Clapping massage	

Note: If you want to use voice control, please say "Hello massage chair" to wake up the chair first, and then use the voice command after getting the reply "Hello Master".

Cleaning and Maintenance


Maintenance methods

Please follow the below steps to clean leather and plastic parts

1. Please ensure the main switch is off and the power plug is unplugged from the electrical outlet.
2. Dip a soft, dry cloth in water or warm water containing 3% to 5% neutral detergent and wring it thoroughly.
3. Clean the leather and plastic parts.
4. Then wipe with dry soft cloth.
5. Let it dry naturally.

Warning


- Do not use alcohol, thinner, gasoline and other chemicals for cleaning.
- Do not dry the surface with a blower.
- Do not iron the fabric.
- During maintenance, always unplug the power plug. Do not unplug it by wet hands.
- Please use a dust cover to avoid dust if you won't use the chair for a long time.
- As fabric's color may be dyed on the leather, be careful when using denim and colored cloth to wipe the product.
- Wipe the product with the melamine foam dipped in mild detergent, if the dirt is hard to remove.
- Do not let the leather touch plastic in case of the color fading.
- Water temperature should not exceed 40 Celsius degrees

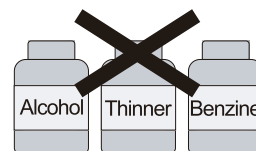Please do not use

Water Temperature < 40°C

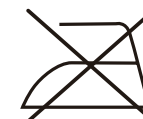Do not iron


Troubleshoot


If there is still an abnormality after the following checks, to prevent any accident, please stop using the massage chair immediately. Please unplug the power cord from the outlet and contact with the dealer or the nearby customer service center.

Problem	Possible cause/action	
The massage rollers move up and then stop.	<ul style="list-style-type: none"> In order not to affect the body, the automatic timer function will move the massage roller to the upper end when the massage is over . 	Inevitable in structure, not abnormal
Noise during the massage	<ul style="list-style-type: none"> The sound of clapping; Friction sound of the massage wheel and fabric during kneading; The sound of the massage wheel when lifting and kneading; The sound during air pressure massage; The sound of the inflator pump; The sound of airbag exhausting; The sound of the backrest up and down. 	
The massage rollers will stop during the massage.	<ul style="list-style-type: none"> If the force applied to the massage rollers is too large, the massage rollers will stop running. Lift your body slightly when using it, or raise the backrest slightly. (especially those who weigh more than 100kg, please pay attention when using) 	After checking, if there is still an abnormality, please stop using the massage chair immediately.
Press the power button, but the massage chair doesn't work properly or cannot connect to the power.	<ul style="list-style-type: none"> Please check the following: <ol style="list-style-type: none"> Is the power plug fully inserted into the socket? Is the main power switch on the left rear of the unit turned on? 	
The power cord and plug are abnormally hot.	<ul style="list-style-type: none"> Please stop using immediately. 	Please contact with your dealer or customer service. Do not disassemble or repair it yourself or by non-professionals.

Specifications

Model Number	Damico Altro 3D
Rated Voltage	100-240V~
Rated Frequency	50-60Hz
Rated Power	170W
Upright position - Approximately	1547mm(L)x750mm(W)x1153mm(H)
Reclined position - Approximately	1700mm(L)x750mm(W)x850mm(H)

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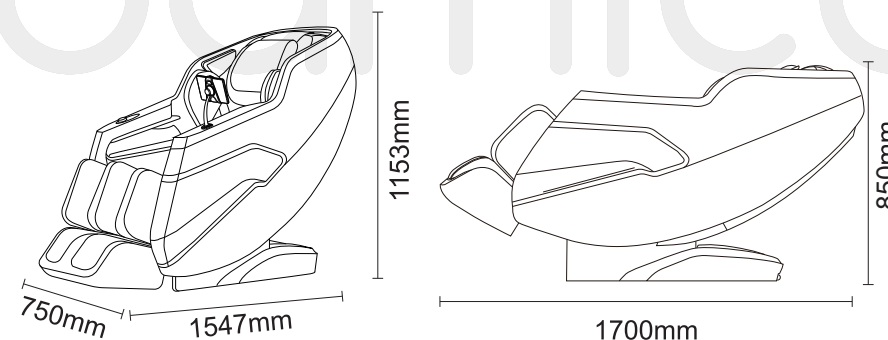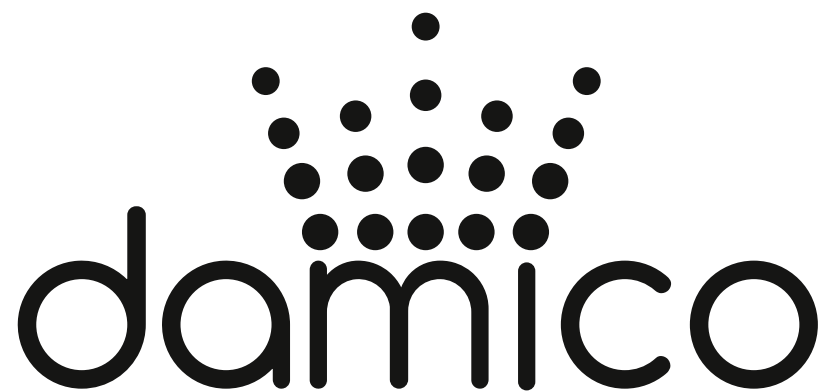Intention	Corpus	Reply	
Basic Functions	Hello massage chair	Hello Master	
	Power on	Please select your favorite program to start your massage.	
	Power off	During shutdown, please leave or restart the massage chair after the chair is completely reset. Thank you and hope you enjoy the massage.	
	Pause massage	Massage paused	
	Continue to massage	Massage resumed	
	Turn on Bluetooth	OK, Bluetooth on	
	Turn off Bluetooth	OK, Bluetooth off	
	Turn off Foot roller	OK, Foot roller off	
	Turn on Foot roller	OK, Foot roller on	
	Turn on Heating	OK, Heating on	
	Turn off Heating	OK, Heating off	
	Zero gravity	OK,Command complete	
	Back up	OK,Command complete	
	Back down	OK,Command complete	
	Leg up	OK,Command complete	
	Leg down	OK,Command complete	
	Start airbag massage	OK, Start airbag massage	
	Stop airbag massage	OK, Stop airbag massage	
	Auto Massage setting	Master massage	Start Master massage
		Relaxing massage	Start Relaxing massage
Deep Shiatsu massage		Start Deep Shiatsu massage	
Refreshing massage		Start Refreshing massage	
Wake Up massage		Start Wake Up massage	
Noon Rest massage		Start Noon Rest massage	
Neck & Shoulder massage		Start Neck & Shoulder massage	
Waist & Back massage		Start Waist & Back massage	
Slim massage		Start Slim massage	

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